

thermo scientific



# Thermo Scientific ClipTip pipetting systems

Sealed in security | Feel the difference

**ThermoFisher**  
SCIENTIFIC

# Breakthrough ClipTip technology

Achieve newfound confidence with the **Thermo Scientific™ ClipTip™ pipetting system** that locks tips firmly in place so they will not loosen or fall off, regardless of application pressure. You can actually feel the tips lock firmly in place with a light touch, making daily pipetting a totally new experience. No more banging tips on your pipette to ensure that they are properly attached and sealed.



**Worry-free pipetting is finally here with the only tip that clips.**

**Eliminate the frustration of banging, checking and rocking tips on your pipette**

With a light touch, each tip is 'clipped' securely on the pipette.

**Minimize the risk of loose or leaking tips**

When ClipTip pipette tips are on, they are sealed and will not loosen. The ClipTip pipetting system remains sealed until you decide to eject tips, even in the most comprehensive applications, including multiple tip touches and mixing cycles.

**Enable consistent and reproducible pipetting from user to user**

Higher-quality results and more efficient research.

4 | **E1-ClipTip Electronic Pipettes**

6 | **F1-ClipTip Manual Pipettes**

12 | **F1-ClipTip System Order Information**

14 | **E1-ClipTip System Order Information**


You'll feel the difference the first time you use the Thermo Scientific ClipTip pipetting system.

## It is a challenge to ensure a quality seal in daily pipetting.

The ideal solution is for the pipette and tip to form a system that increases confidence in reproducibility, reduces forces required to attach and detach the system, and secures the best possible accuracy and precision.



### ➤ FRICTION SEALING SYSTEMS

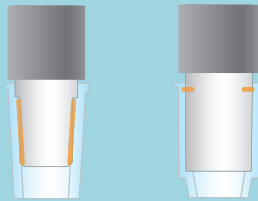



Pipettes with a friction-seal system rely on a user's force to attach tips to the pipette, which varies from user to user.

- Increased tip attachment and ejection forces
- Loose tips that can drop off or leak air
- Decreased confidence in reproducibility of results
- Premature tip cone wear

### ➤ CLIPTIP SYSTEM

By clipping the tips to the pipette with a light touch the ClipTip system optimizes comfort, accuracy and precision

- Extremely low tip attachment and ejection forces
- Reproducibility you can count on
- Complete seal on every channel
- No loose tips

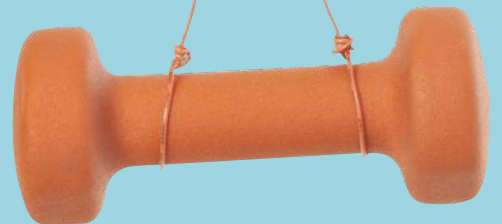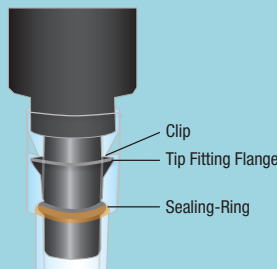

## “They are user-friendly, accurate

and the ‘ClipTip’ technology makes it easy to change tips in the middle of an assay without any difficulties. When using a multichannel its easy for tips to fall off, resulting in inaccurate volumes and if really unlucky, contamination; but with the E1-ClipTip pipettes the tips clip on effortlessly and stay on throughout your sampling. We have various versions of the E1 pipettes in our lab and plan to continue using because of its superior quality!”

– Lab Manager, US

See it in action at:

[thermofisher.com/clipitip](http://thermofisher.com/clipitip)



# Electronic pipettes

Feel the difference with E1-ClipTip Pipettes

### **An optimal solution for complex protocols**

When combined with interlocking ClipTip technology, the electronic pipetting action and flexible user interface of the E1-ClipTip System enable efficient use of precious time and resources.

### **Electronic operation helps you enjoy pipetting again**

Designed for optimum comfort and ease of use, the electronic tip ejection and index finger pipetting action let your thumb relax during pipetting.

### **Ideal operation for multiple users and environments**

Personalized user interface and password protection ensure unified operation among multiple users.

### **E1-ClipTip Equalizer – a perfect choice for sample transfers between various labware formats**

E1-ClipTip Equalizer pipette with adjustable tip spacing allows you to set the distance between tips simply by sliding the scale to expand and contract to your desired setting. This means fewer repetitions for multiple applications.

### **Connected Pipetting via Bluetooth or USB cable**

Connect your E1-ClipTip pipette to the My Pipette Creator app on your PC to easily create programs, to download pre-written protocols, and to share protocols between pipettes and colleagues.

# Enjoy pipetting again

One pipette, one movement, one press of the trigger

## Convenient Li-ion battery

Quickly recharge for long-lasting performance

## Maximized comfort

Electronic index finger pipetting action minimizes thumb strain and an adjustable finger rest allows for optimal pipetting position for both left and right-handed users. Rotating display, unique multichannel ski boot design, and short single-channel barrel offer visibility and a neutral pipetting position in any working environment.

## Program storage

Store up to 30 protocols in your pipette - and the rest in Thermo Fisher Cloud. Programming protocols on your PC with My Pipette Creator application reduces strain on your eyes and thumbs, also allowing you to transfer protocols to multiple pipettes, instead of programming individually.

## Pipette tracking

Service and calibration trackers remind you when the pipette needs to be re-calibrated or sent for service

## Effortless ejection

Electronic tip ejection is effortless, and intuitive user interface prevents any accidental tip ejection whenever there is liquid inside the tips

## A pipette for multiple applications

Choose from single channels, 8–16 channel multichannels, or 6–12 channel adjustable tip spacing multichannel pipettes within the volume range of 0.5–1250  $\mu\text{L}$  based on your application and consumables type

## Customizable calibrations

Higher quality results with customizable calibrations. Improve accuracy by creating different calibrations for liquids with varying viscosity or when using special tips.

## E1-Cliptip pipette

Requires 87% lower tip attachment forces and up to 93% lower tip ejection forces compared to other adjustable tip spacing pipettes<sup>2</sup>



<sup>1</sup> E1-ClipTip electronic pipettes are warranted to be free from defects in material and workmanship for a period of two years from date of purchase. Two-year warranty applies to pipettes that have been registered, one-year warranty to non-registered. To register the pipette and for more information, visit [www.thermoscientific.com/pipettewarranty](http://www.thermoscientific.com/pipettewarranty)

<sup>2</sup> Thermo Scientific Smart Note. Tip attachment and ejection force vary significantly among different adjustable tip spacing pipette brands. 2016.

# Intuitive operation & optimal functions to feel comfortable from day one




## Pipettes that work the way you do

The icon-based user interface offers clear, fast and easy operation with the same basic functions that you are familiar with on your current pipette. My Pipette Function allows you to customize your main menu view to include your most-used functions. You can also create shortcuts for your preferred programs and pipetting functions.


**Now with BLUETOOTH connectivity and 5 user profiles**

## Two programming options:

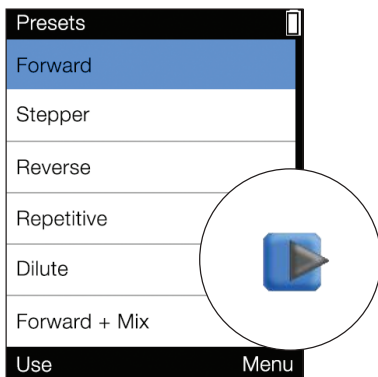


### Presets for everyday ready-to-go mode

Simple and quick pipette settings for your most common applications, such as forward and reverse pipetting and multidispensing. Adjust the volume and speed, and you're ready to go.

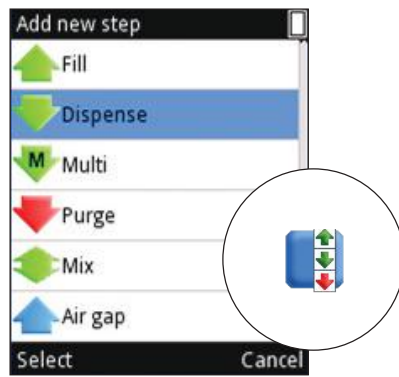

### Matrix for complex applications

The step-based programming is an excellent option when a variety of pipetting functions is required in a unique order. Rather than trying to make a pre-set function work, you can select any of the functional steps — such as Fill, Dispense, Mix, Multi, or Purge — together while setting specific volumes and speeds for each.

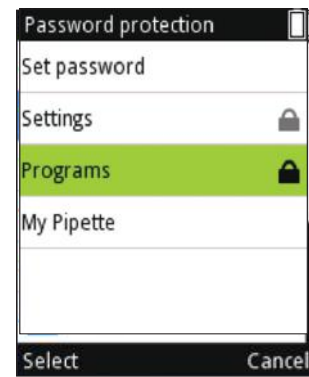

### Choose to lock important protocol settings

Safeguard the integrity of your assays by utilizing password protection for saved settings such as calibration and stored programs, as well as a customized main menu view to deliver consistent, efficient operation regardless of the user.


# Introducing web-based My Pipette™ Creator pipetting app

Boost your productivity and efficiency and connect your E1-ClipTip pipette via Bluetooth™ or USB cable to the My Pipette™ Creator app on your PC. My Pipette™ Creator enables you to easily create programs and download pre-written protocols that can be shared with colleagues and multiple pipettes at once.



## Download protocols

Download pre-programmed protocols for your E1-ClipTip pipettes and favorite Thermo Fisher Scientific reagent kits from the Protocol Library. Save time, enhance consistency, minimize programming errors, and help ensure optimal results.

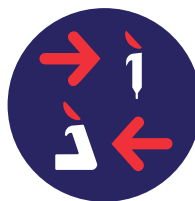

## Program protocols

Reduce the strain on your eyes and thumbs by programming protocols on your PC. Minimize programming errors and enable consistent results by transferring protocols to multiple pipettes, instead of programming each pipette individually.



## Share protocols

Share protocols with your colleagues in another lab, location or country.



## Manage pipettes

Connect multiple E1-ClipTip pipettes to the My Pipette Creator application for efficient, centralized programming that provides consistent results and improved collaboration.

**Learn more at [thermofisher.com/mypipette](https://thermofisher.com/mypipette)**

# Sample transfers have never been this simple

## Do sample transfers require several different pipette models and take up too much time?

Perform sample transfers between virtually any tube, rack, microplate, or horizontal gel box.

Adjustable tip spacing allows you to set the distance between tips simply by sliding the scale to expand and contract to your desired setting.

The unique equalizer window links the tip spacing scale to the particular application. This means fewer repetitions for multiple applications.

## Transfer multiple samples between different lab formats *at once*



### Increased research efficiency and quality





Over **80%** time savings and fewer repetitions compared to single channel pipettes<sup>\*</sup>


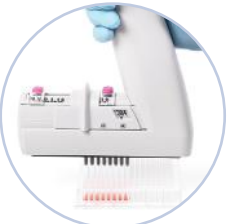




\*Compared to a manual single channel pipette filling a 96-well microplate from microcentrifuge tubes. See Thermo Scientific Smart Note: *Is there an easier and more efficient way to transfer liquids between various labware formats, than using traditional handheld pipettes?* 2014.



# Speed up sample transfer time and reduce repetition

E1-ClipTip Equalizer Pipettes					
Channels	Maximum spacing (mm)	Transfer to/from	↔	Transfer from/to	Minimum spacing (mm)
<b>6-channel</b>	19.8	<ul style="list-style-type: none"> <li>• Test tube racks</li> <li>• 24-well microplates</li> </ul> 	↔	<ul style="list-style-type: none"> <li>• 96-well microplates</li> <li>• Deepwell blocks</li> </ul> 	9.0
<b>8-channel</b>	14.2	<ul style="list-style-type: none"> <li>• Microcentrifuge tube racks</li> <li>• 48-well microplates</li> </ul> 	↔	<ul style="list-style-type: none"> <li>• 96-well microplates</li> <li>• Deepwell blocks</li> </ul> 	9.0

Equalizer Pipettes					
Channels	Maximum spacing (mm)	Transfer to/from	↔	Transfer from/to	Minimum spacing (mm)
<b>8-channel</b>	14.2	<ul style="list-style-type: none"> <li>• 96-well microplates</li> <li>• 48-well microplates</li> <li>• Deepwell blocks</li> <li>• Microcentrifuge tube racks</li> </ul> 	↔	<ul style="list-style-type: none"> <li>• 384-well microplates</li> <li>• 384 PCR plates</li> <li>• Agarose Gels</li> </ul> 	4.5
<b>12-channel</b>	9.0	<ul style="list-style-type: none"> <li>• 96-well microplates</li> <li>• Deepwell blocks</li> </ul> 	↔	<ul style="list-style-type: none"> <li>• 384-well microplates</li> <li>• 384 PCR plates</li> <li>• Agarose Gels</li> </ul> 	4.5

# Manual pipettes

## Maintain precise sample volumes on every channel<sup>1</sup>

Whether you're using a single or multichannel pipette, the F1-ClipTip pipetting system enables a complete seal on every channel for precise sample volumes, saving valuable time to accelerate your discoveries.

## Attach tips easily to help reduce strain injuries<sup>2</sup>

No longer bang tips onto your pipette. With a light touch each tip is 'clipped' securely on the F1-ClipTip pipette and remains sealed until you decide to eject tips.

## A complete tip seal for confidence in your results<sup>1</sup>

The F1-ClipTip system enables consistent and reproducible pipetting from user to user.

## "We love these ClipTip pipettes!"

Hearing the click is an excellent indicator that I have secured the tips to the multichannel. Some of the other brands do not have anything close to that security – many times using the multichannel was impossible because I would have to push each tip on securely, removing any time-savings the pipettor claimed to offer. It's so light weight too!

To date our real-time results have been more consistent than with other multichannel pipettes and so reliably precise we use our set exclusively for real-time now."

– Grace W., Lead Research Specialist



See it in action at: [thermofisher.com/cliptip](https://www.thermofisher.com/cliptip)

<sup>1</sup> Thermo Scientific Application Note: ClipTip Technology Part 2 - Secure Seals for Every Channel, 2012

<sup>2</sup> Thermo Scientific Application Note: ClipTip Technology Part 3 - Reduce the Risk of Repetitive Strain Injury (RSI), 2013

# Exceptional performance and comfort, redefined

## Adjustable comfort

The 120° adjustable finger rest – for both right- and left-handed operator comfort – allows your hand to relax during pipetting cycles.

## Industry-leading lifetime warranty\*

With convenient online registration.

## ClipTip Technology

A secure seal on every channel with low attachment and ejection forces for results you can count on again and again.

## Silent, effortless volume adjustment

Enjoy the new silent and light volume adjustment with over 50% lower forces compared to earlier models. Improved grip on the volume adjustment knob offers effortless and precise adjustment.



## Secure volume lock

Securely set and lock volumes to prevent undesirable volume drift during operation.

\*Thermo Scientific F1-ClipTip pipettes are free from defects in material and workmanship for a product lifetime (up to 10 years). Validity of lifetime warranty varies by country. For more details, visit [www.thermoscientific.com/pipettewarranty](http://www.thermoscientific.com/pipettewarranty)

# F1-ClipTip Manual Pipettes and ClipTip Pipette Tips

Item No.	Description	Color Key	Volume Range	Increment/Packaging
<b>2/10 µL Pipettes</b>				<b>Increment</b>
4641310N	F1-ClipTip 0.1-2 µL	●	0.1-2 µL	0.002 µL
4641320N	F1-ClipTip 1-10 µL	●	1-10 µL	0.02 µL
4661210N	F1-ClipTip 8-ch 1-10 µL	●	1-10 µL	0.02 µL
4661220N	F1-ClipTip 12-ch 1-10 µL	●	1-10 µL	0.02 µL
<b>12.5 Ext µL Tips*</b>				<b>Packaging</b>
94410060	ClipTip 12.5 Ext	●	0.1-12.5 µL	Racked, 96 tips/rack, 10 racks/pack
94410063	ClipTip 12.5 Ext	●	0.1-12.5 µL	Racked, Sterile, 96 tips/rack, 10 racks/pack
94420063	ClipTip 12.5 Ext	●	0.1-12.5 µL	Racked, Filtered, Sterile, 96 tips/rack, 10 racks/pack
<b>20 µL Pipettes</b>				
4641180N	F1-ClipTip 2-20 µL		2-20 µL	0.02 µL
<b>20 µL TIPS †</b>				<b>Packaging</b>
94410210	ClipTip 20	●	1-20 µL	Racked, 96 tips/rack, 10 racks/pack
94410213	ClipTip 20	●	1-20 µL	Racked, Sterile, 96 tips/rack, 10 racks/pack
94410217	ClipTip 20	●	1-20 µL	Reload Tower, 96 tips/insert, 10 inserts/pack
94410218	ClipTip 20	●	1-20 µL	Reload Tower, Sterile, 96 tips/insert, 10 inserts/pack
94420213	ClipTip 20	●	1-20 µL	Racked, Filtered, Sterile, 96 tips/rack, 10 racks/pack
94420218	ClipTip 20	●	1-20 µL	Reload Insert, Filtered, Sterile, 96 tips/insert, 10 inserts/pack
<b>50 µL PIPETTES</b>				<b>Increment</b>
4641190N	F1-ClipTip 5-50 µL	●	5-50 µL	0.1 µL
4661120N	F1-ClipTip 8-ch 5-50 µL	●	5-50 µL	0.1 µL
4661160N	F1-ClipTip 12-ch 5-50 µL	●	5-50 µL	0.1 µL
<b>50 µL TIPS</b>				<b>Packaging</b>
94410250	ClipTip 50	●	5-50 µL	Racked, 96 tips/rack, 10 racks/pack
94410253	ClipTip 50	●	5-50 µL	Racked, Sterile, 96 tips/rack, 10 racks/pack
94410257	ClipTip 50	●	5-50 µL	Reload Tower, 96 tips/insert, 10 inserts/pack
94410258	ClipTip 50	●	5-50 µL	Reload Tower, Sterile, 96 tips/insert, 10 inserts/pack
94420253	ClipTip 50	●	5-50 µL	Racked, Filtered, Sterile, 96 tips/rack, 10 racks/pack
94420258	ClipTip 50	●	5-50 µL	Reload Insert, Filtered, Sterile, 96 tips/insert, 10 inserts/pack
<b>100 µL / 200 µL PIPETTES</b>				<b>Increment</b>
4641200N	F1-ClipTip 10-100 µL	●	10-100 µL	0.2 µL
4641210N	F1-ClipTip 20-200 µL	●	20-200 µL	0.2 µL
4661130N	F1-ClipTip 8-ch 10-100 µL	●	10-100 µL	0.2 µL
4661170N	F1-ClipTip 12-ch 10-100 µL	●	10-100 µL	0.2 µL
<b>200 µL TIPS</b>				<b>Packaging</b>
94410310	ClipTip 200	●	10-200 µL	Racked, 96 tips/rack, 10 racks/pack
94410313	ClipTip 200	●	10-200 µL	Racked, Sterile, 96 tips/rack, 10 racks/pack
94410317	ClipTip 200	●	10-200 µL	Reload Tower, 96 tips/insert, 10 inserts/pack
94410318	ClipTip 200	●	10-200 µL	Reload Tower, Sterile, 96 tips/insert, 10 inserts/pack
94420313	ClipTip 200	●	10-200 µL	Racked, Filtered, Sterile, 96 tips/rack, 10 racks/pack
94420318	ClipTip 200	●	10-200 µL	Reload Insert, Filtered, Sterile, 96 tips/insert, 10 inserts/pack
<b>300 µL PIPETTES</b>				<b>Increment</b>
4641220N	F1-ClipTip 30-300 µL	●	30-300 µL	1 µL
4661140N	F1-ClipTip 8-ch 30-300 µL	●	30-300 µL	1 µL
4661180N	F1-ClipTip 12-ch 30-300 µL	●	30-300 µL	1 µL
<b>300 µL TIPS</b>				<b>Packaging</b>
94410510	ClipTip 300	●	30-300 µL	Racked, 96 tips/rack, 10 racks/pack
94410513	ClipTip 300	●	30-300 µL	Racked, Sterile, 96 tips/rack, 10 racks/pack
94410517	ClipTip 300	●	30-300 µL	Reload Tower, 96 tips/insert, 10 inserts/pack

- 2 µL/10 µL/12.5 µL/20 µL
- 50 µL
- 100 µL/ 200 µL
- 300 µL
- 850 µL/ 1000 µL
- 1250 µL

\* CT 12.5 Ext tips and compatible F1-ClipTip pipettes utilize 'snap and seal' technology. All ClipTip pipette tips are low retention except for the CT12.5, CT12.5 Ext and 384-format tips.  
 † CT20 tips are also used with old F1-ClipTip 1-10µl pipette models with catalogue numbers 4641170N, 4661110N, 4661150N and 4651190N.

94410518	ClipTip 300	●	30-300 µL	Reload Tower, Sterile, 96 tips/insert, 10 inserts/pack
94420513	ClipTip 300	●	30-300 µL	Racked, Filtered, Sterile, 96 tips/rack, 10 racks/pack
94420518	ClipTip 300	●	30-300 µL	Reload Insert, Filtered, Sterile, 96 tips/inset, 10 inserts/pack
94410610	ClipTip 300 Ext	●	30-300 µL	Racked, 96 tips/rack, 8 racks/pack
94410613	ClipTip 300 Ext	●	30-300 µL	Racked, sterile, 96 tips/rack, 8 racks/pack
94410617	ClipTip 300 Ext	●	30-300 µL	Reload tower, 96 tips/insert, 10 inserts/pack
94410618	ClipTip 300 Ext	●	30-300 µL	Reload tower, sterile, 96 tips/insert, 10 inserts/pack
94420613	ClipTip 300 Ext	●	30-300 µL	Racked, filtered, sterile, 96 tips/rack, 8 racks/pack

● 2 µL/10 µL/12.5 µL/20 µL  
 ● 50 µL  
 ● 100 µL/ 200 µL  
 ● 300 µL  
 ● 1000 µL  
 ● 1250 µL

Item No.	Description	Color Key	Volume Range	Increment/Packaging
<b>1000 µL PIPETTES</b>				<b>Increment</b>
4641230N	F1-ClipTip 100-1000 µl**	●	100-1000 µL	1 µL
<b>1000 µL TIPS</b>				<b>Packaging</b>
94410710	ClipTip 1000	●	100-1000 µL	Racked, 96 tips/rack, 8 racks/pack
94410713	ClipTip 1000	●	100-1000 µL	Racked, Sterile, 96 tips/rack, 8 racks/pack
94410717	ClipTip 1000	●	100-1000 µL	Reload Tower, 96 tips/insert, 8 inserts/pack
94410718	ClipTip 1000	●	100-1000 µL	Reload Tower, Sterile, 96 tips/insert, 8 inserts/pack
94420713	ClipTip 1000	●	100-1000 µL	Racked, Filtered, Sterile, 96 tips/rack, 8 racks/pack

Item No.	Description	Volume Range
4651280N	F1-ClipTip 10 µL	10 µL
4651200N	F1-ClipTip 20 µL	20 µL
4651210N	F1-ClipTip 25 µL	25 µL
4651220N	F1-ClipTip 50 µL	50 µL
4651230N	F1-ClipTip 100 µL	100 µL
4651240N	F1-ClipTip 200 µL	200 µL
4651250N	F1-ClipTip 250 µL	250 µL
4651260N	F1-ClipTip 500 µl**	500 µL
4651270N	F1-ClipTip 1000 µl**	1000 µL

Order online and find additional specifications at [thermofisher.com/cliptip](https://www.thermofisher.com/cliptip)

\*\* Short Shaft

Item No.	Description	Contents of Kit
4701140N	F1-ClipTip GLP Kit 1	<ul style="list-style-type: none"> <li>• F1-ClipTip 1-10 µL single channel</li> <li>• F1-ClipTip 10-100 µL single channel</li> <li>• F1-ClipTip 100-1000 µL single channel</li> <li>• F-Series Stand</li> <li>• ClipTip 20 rack</li> <li>• ClipTip 200 rack</li> <li>• ClipTip 1000 rack</li> <li>• 25 mL reservoir</li> <li>• 25 mL divided reservoir</li> <li>• 100 mL divided reservoir</li> </ul>
4701150N	F1-ClipTip GLP Kit 2	<ul style="list-style-type: none"> <li>• F1-ClipTip 2-20 µL single channel</li> <li>• F1-ClipTip 20-200 µL single channel</li> <li>• F1-ClipTip 100-1000 µL single channel</li> <li>• F1-ClipTip 30-300 µL 8-channel</li> <li>• F-Series Stand</li> <li>• ClipTip 20 rack</li> <li>• ClipTip 200 rack</li> <li>• ClipTip 300 rack</li> <li>• ClipTip 1000 rack</li> <li>• 25 mL reservoir</li> <li>• 25 mL divided reservoir</li> <li>• 100 mL divided reservoir</li> </ul>



Item No.	Description	Color Key	Increment/Packaging
9420400	F-Series Stand, white, 6 positions		
94410219	Empty Rack, Small	● ●	10 racks, no tips/inserts
94410519	Empty Rack, Medium	● ●	10 racks, no tips/inserts
94410619	300 Ext Empty Rack	●	8 racks, no tips/inserts
94410819	Empty Rack, Large	● ●	8 racks, no tips/inserts

# E1-ClipTip Electronic Pipettes and ClipTip Pipette Tips

- 2 µL/10 µL/12.5 µL/20 µL
- 30 µL/ 50 µL
- 100 µL/ 125 µL/ 200 µL
- 300 µL
- 850 µL/ 1000 µL
- 1250 µL

Item No.	Description	Color Key	Volume Range	Tip Spacing Range (mm)	Increment/Packaging
<b>0.5-12.5 µL PIPETTES*</b>					<b>Increment</b>
4670000BT	E1-ClipTip, 0.5-12.5 µL	<span style="color: #e91e63;">●</span>	0.5-12.5 µL		0.01 µL
4671000BT	E1-ClipTip, 8-ch 0.5-12.5 µL	<span style="color: #e91e63;">●</span>	0.5-12.5 µL		0.01 µL
4671010BT	E1-ClipTip, 12-ch 0.5-12.5 µL	<span style="color: #e91e63;">●</span>	0.5-12.5 µL		0.01 µL
<b>12.5 µL TIPS*</b>					<b>Packaging</b>
94410040	ClipTip 12.5	<span style="color: #e91e63;">●</span>	0.5-12.5 µL		Racked, 96 tips/rack, 10 racks/pack
94410043	ClipTip 12.5	<span style="color: #e91e63;">●</span>	0.5-12.5 µL		Racked, sterile, 96 tips/rack, 10 racks/pack
94420043	ClipTip 12.5	<span style="color: #e91e63;">●</span>	0.5-12.5 µL		Racked, filtered, sterile, 96 tips/rack, 10 racks/pack
94410060	ClipTip 12.5 Ext	<span style="color: #e91e63;">●</span>	0.5-12.5 µL		Racked, 96 tips/rack, 10 racks/pack
94410063	ClipTip 12.5 Ext	<span style="color: #e91e63;">●</span>	0.5-12.5 µL		Racked, sterile, 96 tips/rack, 10 racks/pack
94420063	ClipTip 12.5 Ext	<span style="color: #e91e63;">●</span>	0.5-12.5 µL		Racked, filtered, sterile, 96 tips/rack, 10 racks/pack
<b>384, 12.5 µL PIPETTES*</b>					<b>Increment</b>
4671020BT	E1-ClipTip, 16-ch 0.5-12.5 µL	<span style="color: #e91e63;">●</span>	0.5-12.5 µL		0.01 µL
4672010BT	E1-ClipTip Equalizer 384, 8-ch 0.5-12.5 µL	<span style="color: #e91e63;">●</span> Eq	0.5-12.5 µL	4.5-14.2	0.01 µL
4672020BT	E1-ClipTip Equalizer 384, 12-ch 0.5-12.5 µL	<span style="color: #e91e63;">●</span> Eq	0.5-12.5 µL	4.5-9.0	0.01 µL
<b>384, 12.5 µL TIPS*</b>					<b>Packaging</b>
94410050	ClipTip 384, 12.5	<span style="color: #e91e63;">●</span>	0.5-12.5 µL		Racked, 384 tips/rack, 10 racks/pack
94410053	ClipTip 384, 12.5	<span style="color: #e91e63;">●</span>	0.5-12.5 µL		Racked, sterile, 384 tips/rack, 10 racks/pack
94420053	ClipTip 384, 12.5	<span style="color: #e91e63;">●</span>	0.5-12.5 µL		Racked, filtered, sterile, 384 tips/rack, 10 racks/pack
94410070	ClipTip 384, 12.5 Ext	<span style="color: #e91e63;">●</span>	0.5-12.5 µL		Racked, 384 tips/rack, 10 racks/pack
94410073	ClipTip 384, 12.5 Ext	<span style="color: #e91e63;">●</span>	0.5-12.5 µL		Racked, sterile, 384 tips/rack, 10 racks/pack
94420073	ClipTip 384, 12.5 Ext	<span style="color: #e91e63;">●</span>	0.5-12.5 µL		Racked, filtered, sterile, 384 tips/rack, 10 racks/pack
<b>384, 30 µL PIPETTES*</b>					<b>Increment</b>
4671030BT	E1-ClipTip, 16-ch 1-30 µL	<span style="color: #3f51b5;">●</span>	1-30 µL		0.01 µL
4672030BT	E1-ClipTip Equalizer 384, 8-ch 1-30 µL	<span style="color: #3f51b5;">●</span> Eq	1-30 µL	4.5-14.2	0.01 µL
4672040BT	E1-ClipTip Equalizer 384, 12-ch 1-30 µL	<span style="color: #3f51b5;">●</span> Eq	1-30 µL	4.5-9.0	0.01 µL
<b>384, 30 µL TIPS*</b>					<b>Packaging</b>
94410100	ClipTip 384, 30	<span style="color: #3f51b5;">●</span>	1-30 µL		Racked, 384 tips/rack, 10 racks/pack
94410103	ClipTip 384, 30	<span style="color: #3f51b5;">●</span>	1-30 µL		Racked, Sterile, 384 tips/rack, 10 racks/pack
94420103	ClipTip 384, 30	<span style="color: #3f51b5;">●</span>	1-30 µL		Racked, Filtered, Sterile, 384 tips/rack, 10 racks/pack
<b>125 µL PIPETTES</b>					<b>Increment</b>
4670020BT	E1-ClipTip, 2-125 µL	<span style="color: #ffc107;">●</span>	2-125 µL		0.1 µL
4671040BT	E1-ClipTip, 8-ch 2-125 µL	<span style="color: #ffc107;">●</span>	2-125 µL		0.1 µL
4671050BT	E1-ClipTip, 12-ch 2-125 µL	<span style="color: #ffc107;">●</span>	2-125 µL		0.1 µL
4672050BT	E1-ClipTip Equalizer, 8-ch 2-125 µL	<span style="color: #ffc107;">●</span> Eq	2-125 µL	9.0-14.2	0.1 µL
<b>200 µL TIPS</b>					<b>Packaging</b>
94410310	ClipTip 200	<span style="color: #ffc107;">●</span>	2-200 µL		Racked, 96 tips/rack, 10 racks/pack
94410313	ClipTip 200	<span style="color: #ffc107;">●</span>	2-200 µL		Racked, Sterile, 96 tips/rack, 10 racks/pack
94410317	ClipTip 200	<span style="color: #ffc107;">●</span>	2-200 µL		Reload Tower, 96 tips/insert, 10 inserts/pack
94410318	ClipTip 200	<span style="color: #ffc107;">●</span>	2-200 µL		Reload Tower, Sterile, 96 tips/insert, 10 inserts/pack
94420313	ClipTip 200	<span style="color: #ffc107;">●</span>	2-200 µL		Racked, Filtered, Sterile, 96 tips/rack, 10 racks/pack
94420318	ClipTip 200	<span style="color: #ffc107;">●</span>	2-200 µL		Reload Insert, Filtered, Sterile, 96 tips/insert, 10 inserts/pack
<b>384, 125 µL PIPETTES*</b>					<b>Increment</b>
4671060BT	E1-ClipTip, 16-ch 2-125 µL	<span style="color: #ffc107;">●</span>	2-125 µL		0.1 µL
4672060BT	E1-ClipTip Equalizer 384, 8-ch 2-125 µL	<span style="color: #ffc107;">●</span> Eq	2-125 µL	4.5-14.2	0.1 µL
4672070BT	E1-ClipTip Equalizer 384, 12-ch 2-125 µL	<span style="color: #ffc107;">●</span> Eq	2-125 µL	4.5-9.0	0.1 µL
<b>384, 125 µL TIPS*</b>					<b>Packaging</b>
94410150	ClipTip 384, 125	<span style="color: #ffc107;">●</span>	2-125 µL		Racked, 384 tips/rack, 10 racks/pack

94410153	ClipTip 384, 125	●	2-125 µL		Racked, Sterile, 384 tips/rack, 10 racks/pack
94420153	ClipTip 384, 125	●	2-125 µL		Racked, Filtered, Sterile, 384 tips/rack, 10 racks/pack
Item No.	Description	Color Key	Volume Range	Tip Spacing Range (mm)	Increment/Packaging
<b>300 µL PIPETTES</b>					<b>Increment</b>
4670030BT	E1-ClipTip, 10-300 µL	●	10-300 µL		0.1 µL
4671070BT	E1-ClipTip, 8-ch 10-300 µL	●	10-300 µL		0.1 µL
4671080BT	E1-ClipTip, 12-ch 10-300 µL	●	10-300 µL		0.1 µL
4672080BT	E1-ClipTip Equalizer, 8-ch 10-300 µL	Eq ●	10-300 µL	9.0-14.2	0.1 µL
<b>300 µL TIPS</b>					<b>Packaging</b>
94410510	ClipTip 300	●	30-300 µL		Racked, 96 tips/rack, 10 racks/pack
94410513	ClipTip 300	●	30-300 µL		Racked, sterile, 96 tips/rack, 10 racks/pack
94410517	ClipTip 300	●	30-300 µL		Reload tower, 96 tips/insert, 10 inserts/pack
94410518	ClipTip 300	●	30-300 µL		Reload tower, sterile, 96 tips/insert, 10 inserts/pack
94420513	ClipTip 300	●	30-300 µL		Racked, filtered, sterile, 96 tips/rack, 10 racks/pack
94420518	ClipTip 300	●	30-300 µL		Racked, filtered, sterile, 96 tips/rack, 10 racks/pack
94410610	ClipTip 300 Ext	●	30-300 µL		Racked, 96 tips/rack, 8 racks/pack
94410613	ClipTip 300 Ext	●	30-300 µL		Racked, sterile, 96 tips/rack, 8 racks/pack
94410617	ClipTip 300 Ext	●	30-300 µL		Reload tower, 96 tips/insert, 10 inserts/pack
94410618	ClipTip 300 Ext	●	30-300 µL		Reload tower, sterile, 96 tips/insert, 10 inserts/pack
94420613	ClipTip 300 Ext	●	30-300 µL		Racked, filtered, sterile, 96 tips/rack, 8 racks/pack
<b>850 µL PIPETTES</b>					<b>Increment</b>
4671090BT	E1-ClipTip, 12-ch 30-850 µL	●	30-850 µL		1.0 µL
<b>1000 µL TIPS</b>					<b>Packaging</b>
94410710	ClipTip 1000	●	30-1000 µL		Racked, 96 tips/rack, 8 racks/pack
94410713	ClipTip 1000	●	30-1000 µL		Racked, Sterile, 96 tips/rack, 8 racks/pack
94410717	ClipTip 1000	●	30-1000 µL		Reload Tower, 96 tips/insert, 8 inserts/pack
94410718	ClipTip 1000	●	30-1000 µL		Reload Tower, Sterile, 96 tips/insert, 8 inserts/pack
94420713	ClipTip 1000	●	30-1000 µL		Racked, Filtered, Sterile, 96 tips/rack, 8 racks/pack
94420718	ClipTip 1000	●	30-1000 µL		Reload Insert, Filtered, Sterile, 96 tips/insert, 8 inserts/pack
<b>15-1250 µL PIPETTES</b>					<b>Increment</b>
4670040BT	E1-ClipTip, 15-1250 µL	●	15-1250 µL		1.0 µL
4671100BT	E1-ClipTip, 8-ch 15-1250 µL	●	15-1250 µL		1.0 µL
4672090BT	E1-ClipTip Equalizer, 6-ch 15-1250 µL	Eq ●	15-1250 µL	9.0-19.8	1.0 µL
4672100BT	E1-ClipTip Equalizer, 8-ch 15-1250 µL	Eq ●	15-1250 µL	9.0-14.2	1.0 µL
<b>1250 µL TIPS</b>					<b>Packaging</b>
94410810	ClipTip 1250	●	15-1250 µL		Racked, 96 tips/rack, 8 racks/pack
94410813	ClipTip 1250	●	15-1250 µL		Racked, Sterile, 96 tips/rack, 8 racks/pack
94410817	ClipTip 1250	●	15-1250 µL		Reload Insert, 96 tips/insert, 8 inserts/pack
94410818	ClipTip 1250	●	15-1250 µL		Reload Insert, Sterile, 96 tips/insert, 8 inserts/pack
94420813	ClipTip 1250	●	15-1250 µL		Racked, Filtered, Sterile, 96 tips/rack, 8 racks/pack
94420818	ClipTip 1250	●	15-1250 µL		Reload Insert, Filtered, Sterile, 96 tips/insert, 8 inserts/pack

- 2 µL/10 µL/12.5 µL/20 µL
- 30 µL/ 50 µL
- 100 µL/ 125 µL/ 200 µL
- 300 µL
- 850 µL/ 1000 µL
- 1250 µL

Item No.	Description	Color Key	Increment/Packaging
9420500	E1-ClipTip Charging Stand, 1 position		Compatible with E1-ClipTip pipettes
9420510	E1-ClipTip Multi Charging Stand, 3 position		Compatible with E1-ClipTip pipettes
94410219	Empty Rack, Small	● ●	10 racks, no tips/inserts
94410519	Empty Rack, Medium	● ●	10 racks, no tips/inserts
94410619	300 Ext Empty Rack	●	8 racks, no tips/inserts
94410819	Empty Rack, Large	● ●	8 racks, no tips/inserts

Order online and find additional specifications [thermofisher.com/cliptip](http://thermofisher.com/cliptip)

\* ClipTip 384 and ClipTip 12.5 tips and compatible E1-ClipTip models do not have interlock technology. All ClipTip pipette tips are low retention with the exception of 384 and 12.5 models.



E1-ClipTip  
Single  
Charging  
Stand



E1-ClipTip  
Multi  
Charging  
Stand

## Pipetting for a lifetime

A history of liquid handling innovations providing consistent, reproducible results



TECHNIQUES



TOOLS



ERGONOMICS



SERVICE

Learn how the right combination of technique, selecting the right tools, ergonomics and service support Good Laboratory Pipetting at [thermofisher.com/pipettingforlifetime](https://www.thermofisher.com/pipettingforlifetime)

## Protecting your investment

For more than 25 years, we have prided ourselves by offered quality-minded service that meets customers diverse needs.

Whether you're working in a unique environment such as cGLP, cGMP, critical research, academia, or have a special requirement, you can feel confident that your pipettes will have the accuracy and performance you need.

We provide you with:

- Online visibility of service status, calibration history, and asset management
- Tailored solutions to fit your unique needs
- Option to send pipettes to our calibration laboratory or have us come to your site



- **Support Plans**
- **Parts and Consumables**
- **Training and More**



### Let our expertise complement yours

Explore all the new ways to protect your instruments and equipment at [thermofisher.com/equipment-service](https://www.thermofisher.com/equipment-service)

## Pipette warranty registration

Register your pipette online to activate your free warranty at [www.thermofisher.com/pipetwarranty](https://www.thermofisher.com/pipetwarranty)

Find out more at [thermofisher.com/cliptip](https://www.thermofisher.com/cliptip)

© 2021 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific and its subsidiaries unless otherwise specified. The Bluetooth® is a trademark of Bluetooth SIG, Inc. Specifications, terms and pricing are subject to change. Validity of lifetime warranty varies by country. Not all products are available in all countries. Please consult your local sales representative for details.

**ThermoFisher**  
SCIENTIFIC