

smarter FT-IR definitive answers



Nicolet iS50 FT-IR Spectrometer materials analysis

beyond imagination



POLYMERS and RUBBERS

- run one-touch raw materials ID, blending, morphology, and kinetics
- multi-task in a single workstation one setup, multiple answers



FORENSICS

- access two Class-A SWGDRUG methods in a small-footprint unit - simplify GC-IR analysis of isomeric materials to one operation



PHARMACEUTICALS

- analyze formulations and polymorphs via an integrated system
- validate and deploy QA/QC methods for process or loading-dock use

PIGMENTS, PAINTS and INKS

- find inorganic fillers and additives in an organic base
- correlate visible colors with formulation

FOOD, FLAVORS and OILS

- identify isomeric fragrances and flavors - analyze moisture, fat, proteins and extracted oils from seeds

GENERAL LABORATORIES

- bring one-touch consistency to complex analyses
- increase laboratory power with many techniques in a small footprint

The Thermo Scientific Nicolet iS50 FT-IR spectrometer gives definitive answers in your analytical laboratory where fast, easy workflows are essential. No other FT-IR goes beyond your expectations to offer one-touch simplicity, flexibility and integration in such a small, cost-effective workstation. The Nicolet" iS"50 is so advanced you won't think of it as simply an FT-IR.

Go beyond FT-IR to answers.












beyond the ordinary

raise productivity while reducing manual operations

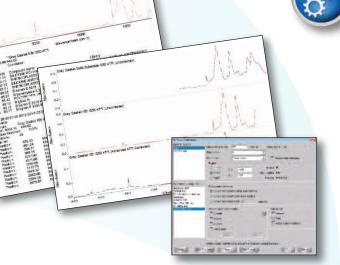

Touch Point operation: collect, analyze and report – you do science, not data processing!



Collect, analyze and compile reports with one touch simplicity

As you add more capability, the Touch Points become more powerful, driving setup and checking performance so you don't need to.

- Select beam path, optical components and sampling module with one touch
- Create automated workflows with no programming
- Compile multiple spectra, analyses and images into one group file
- Track, organize and output information with auto-reporting



Engage your most powerful sampling technique anytime with built-in ATR


ATR like it was meant to be. Rugged, fast and always ready to go. Powders, liquids, gels and pastes can be analyzed, with no setup, using the iS50 ATR. All reflective optics and proprietary detector allow you to analyze organic and inorganic materials with mid- and far-IR spectroscopy. The iS50 ATR leaves the main sample compartment free for other analyses, never requiring you to swap accessories or break system seals.

- Eliminate manual set-up procedures and instrument conditioning time
- Keep the main sample compartment clear for other experiments, such as Raman
- Examine hard, acidic or volatile materials with minimal sample preparation
- Analyze bulk polymers and fillers simultaneously

Automated Multi-range

Redefine analytical productivity with automatic beamsplitter changing


Automate your analysis with the iS50 ABX Automated Beamsplitter exchanger, which enables the Nicolet iS50 to set up the instrument with a single press of a Touch Point. Now, you can do the complete analysis including crystallinity, bulk analysis, and identification, using near-, mid- and far-IR by loading samples and letting the instrument do all the work.

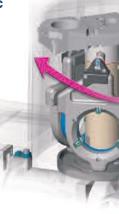

- Measure multiple aspects of your sample automatically
- Enable setup changes without disturbing seals, desiccation or purge
- Eliminate manual handling of optical components



IS50 FT-IR

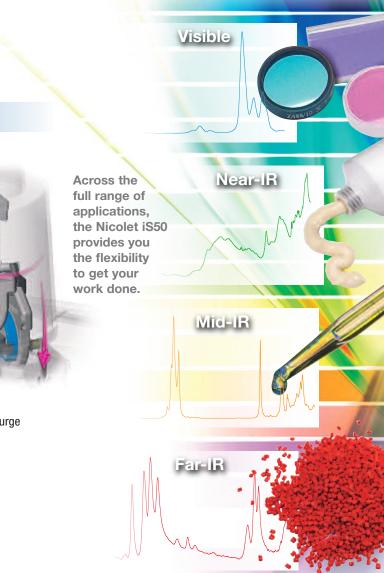

Smart Flexibility

Target your analysis with the right accessory

Smart means getting the right tool for the job. It also means being able to use accessories you already have, with a spectrometer made to accept them.

- Get even more from your Smart Accessories through the Touch Points
- Use a full range of standard accessories
- Full-size sample compartment for large accessories

beyond FT-IR expand your analysis horizons

Integrated Raman


SSOFT-IR

Add the power of Raman to complement your material characterization

The iS50 Raman accessory provides capabilities, like mapping and well-plate screening in a sample compartment mounted module. Never before has laboratory Raman analysis been so versatile, cost effective and easy to use.

- Perform analysis of crystallinity and inorganics
- Complement your material identification with Raman
- Switch in seconds between the iS50 ATR and iS50 Raman
- Screen multiple samples for cluster analysis



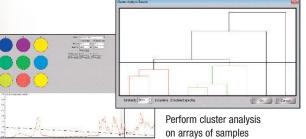



Dedicated NIR

Develop QC methods for bulk samples

The iS50 NIR module provides integrating sphere and fiber optic options in one compact unit.

- Initiate collection with one touch
- Run samples through containers without preparation
- Obtain consistent spectra from heterogeneous solids
- Sample solids and liquids in situ






(

Thermo

Examine powders, tablets, films, fibers or particles

ISSO RAMA




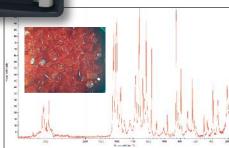

Complete Analysis in One Step

The Nicolet iS50 lets you automate your entire process even when combining techniques. Collect and analyze Raman, mid-IR ATR, near-IR and more in one Macros\Basic driven process. Walk away and come back to a combined report.

- Store images along with your data file
- Couple multiple data files and results for a project
- Build workflows for consistent results with any operator

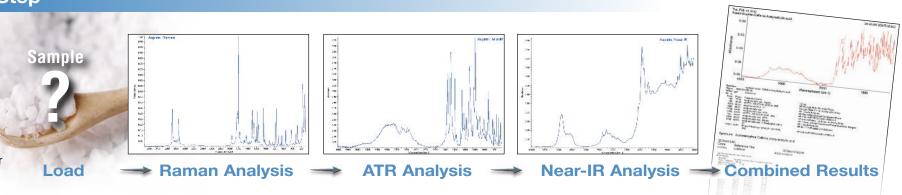



beyond results


bring simplicity and consistency to your operations

Definitive TGA-IR

Make deformulation a routine procedure for your materials laboratory

Complex gas mixtures evolved from the TGA experiment are automatically resolved and identified by the Mercury TGA software.

- Extract "what" and "how much" information from your material
- Get fast, consistent results easily with Mercury TGA software
- Probe starting and residual material via ATR to complete the picture

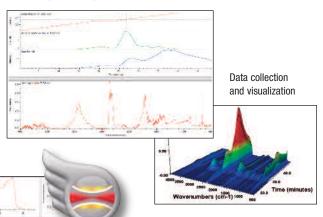

Mercury TGA means TGA-IR is now a genuine analytical tool for any analyst in any laboratory. Put TGA-IR to work.



Thermo

Compare outputs: the differences tell the story for failure analysis or reverse engineering

250 FFH





Streamlined GC-IR


Uniquely identify stereochemistry and trace chemicals

Get qualitative, isomer-specific analysis with no effort using the iS50 GC-IR modules and the Thermo Scientific OMNIC Mercury GC software. Your productivity can rise further using the industry standard Thermo Scientific Dionex Chromeleon chromatography data system software to drive autosamplers and simultaneous GC-FID and GC-IR detection. Add other features of the Nicolet iS50 to create the ultimate forensics laboratory workstation -Raman, ATR and GC - all integrated with the spectrometer.

- Analyze fuels, chemicals and drugs definitively
- Present data meaningfully using waterfall or 3-dimensional plots
- · Extract peaks and identify them automatically using Mercury GC software



As your needs grow, so can the capabilities of the Nicolet iS50 spectrometer. Whether you are studying microscopic structures, organic thin films, fast reactions or new light sources, the Nicolet iS50 spectrometer enables your research now and in the future.

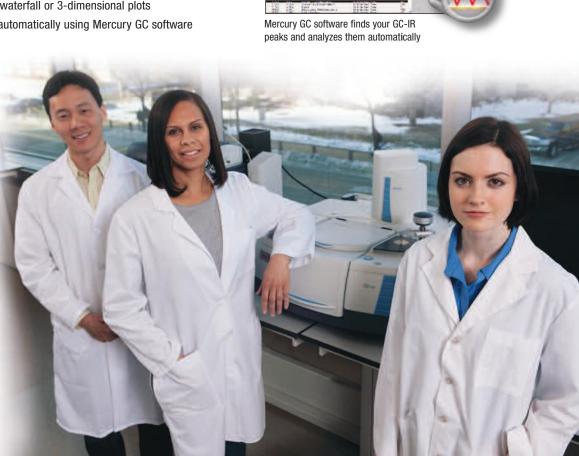

Isomeric materials display unique spectra for identification

beyond tomorrow

the most trusted FT-IR in the industry is...



backed by world-class service and applications expertise. Equipped with tools to facilitate regulatory compliance. Field upgradable to ensure best fit to your future needs.

Our commitment to you ensures your ability to provide definitive answers...today, tomorrow and beyond.



New Zealand +64 9 980 6700 Russia/CIS +43 1 333 50 34 0 Spain +34 914 845 965 Switzerland +41 61 716 77 00 UK +44 1442 233555 USA +1 800 532 4752



thermoscientific.com

©2012 Thermo Fisher Scientific Inc. All rights reserved. ISO is a trademark of the International Standards Organization. All other trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details.

 Africa
 +27
 11
 822
 4120
 Denmark
 +45
 70
 23
 62
 60
 India
 +91
 22
 6742
 94

 Australia
 +61
 3
 9757
 4300
 Europe-Other
 +43
 1
 333
 50
 40
 Italy
 +39
 02
 950
 591
Austria +43 1 333 50 34 0 Finland/Norway/Sweden Belgium +32 53 73 42 41 Canada +1 800 530 8447 **China** +86 10 8419 3588

BR52315_E 04/12M

+46 8 556 468 00 France +33 1 60 92 48 00 Germany +49 6103 408 1014

India +91 22 6742 9434 Japan +81 45 453 9100 Latin America +1 561 688 8700 Middle East +43 1 333 50 34 0 **Netherlands** +31 76 579 55 55